



CHROMATOGRAPHY SIMPLIFIED

When accuracy and precision are everything

The **PerkinElmer LC 300 Platform**, with its wide range of detectors, accessories and the SimplicityChrom CDS software, redefines your liquid chromatography workflow, offering the flexibility, performance and efficiency needed to tackle even the most challenging analytical demands.

For more information about the
LC 300 systems, click here:

LEARN MORE



JOURNAL OF SEPARATION SCIENCE




3 | 2022

Jigucao capsule
A famous Chinese patent medicine



LC-MS/MS

Qualitative Analysis

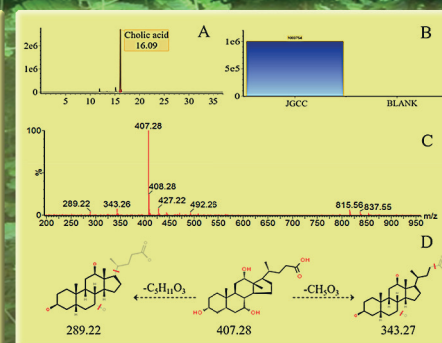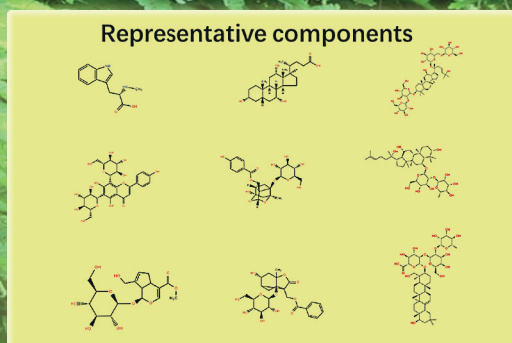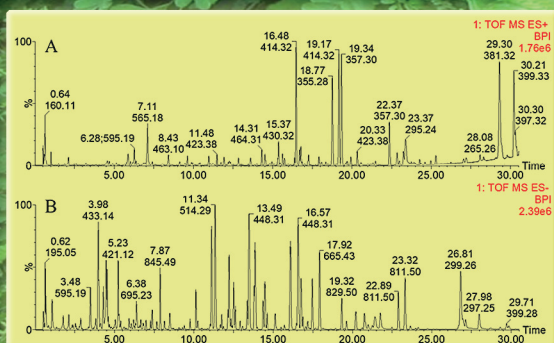

Methods

Chromatography · Electroseparation

Applications

Biomedicine · Foods · Environment

www.jss-journal.com

WILEY-VCH

CORRIGENDUM

Advanced data preprocessing for Comprehensive two-dimensional Gas Chromatography with Vacuum Ultraviolet Spectroscopy detection

Unfortunately, in the article

Advanced data preprocessing for Comprehensive two-dimensional Gas Chromatography with Vacuum Ultraviolet Spectroscopy detection

By Aleksandra Lelevic, Vincent Souchon, Christophe Geantet, Chantal Lorentz, and Maxime Moreaud

J. Sep. Sci. 2021, 44, 4141–4150

DOI: 10.1002/jssc.202100528

Equations (3), (5) and (6) were given incorrectly. The correct Equations (3), (5) and (6) are as follows:

$$F_{\alpha}(s) = s * f_{\alpha} - \text{ReLU}(s * f_{\alpha} - s) \quad (3)$$

$$b = (F_{\alpha} \circ F_{\alpha} \circ \dots \circ F_{\alpha})(s) * f_{\alpha} \quad (5)$$

$$b = (F_{\alpha})^n(s) * f_{\alpha} \quad (6)$$

The changes to these equations lead to slight modifications of Figures 2 and 3. The correct Figures 2 and 3 are as follows:

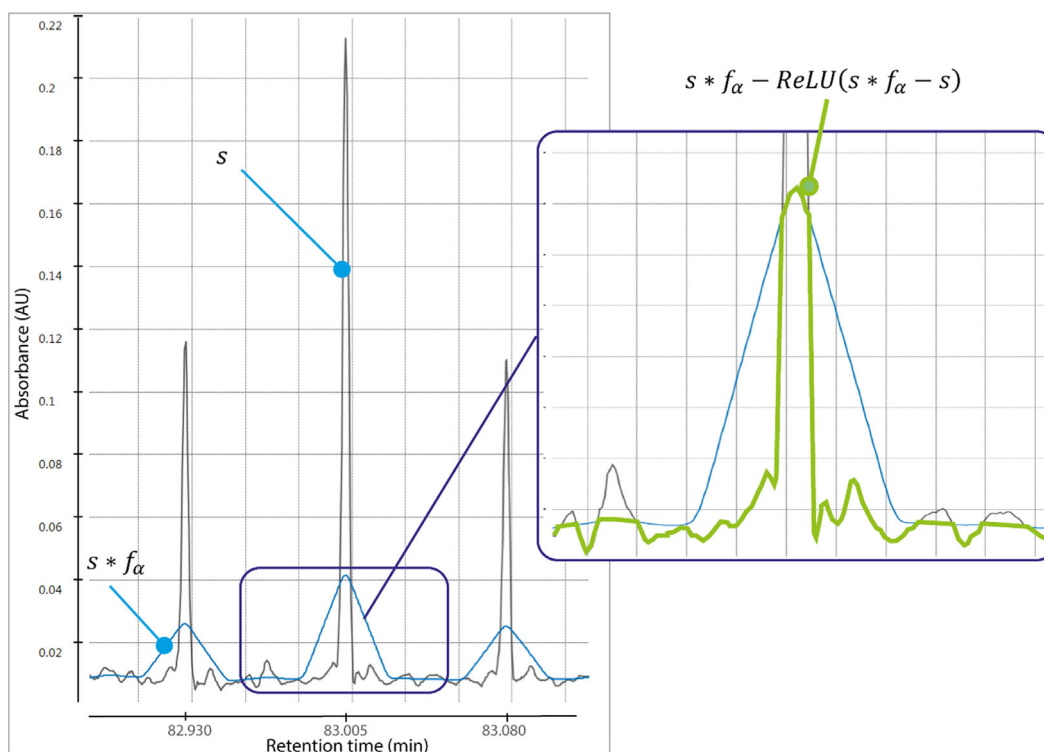


FIGURE 2 Illustration of the approach for the baseline correction ($\alpha = 20$)

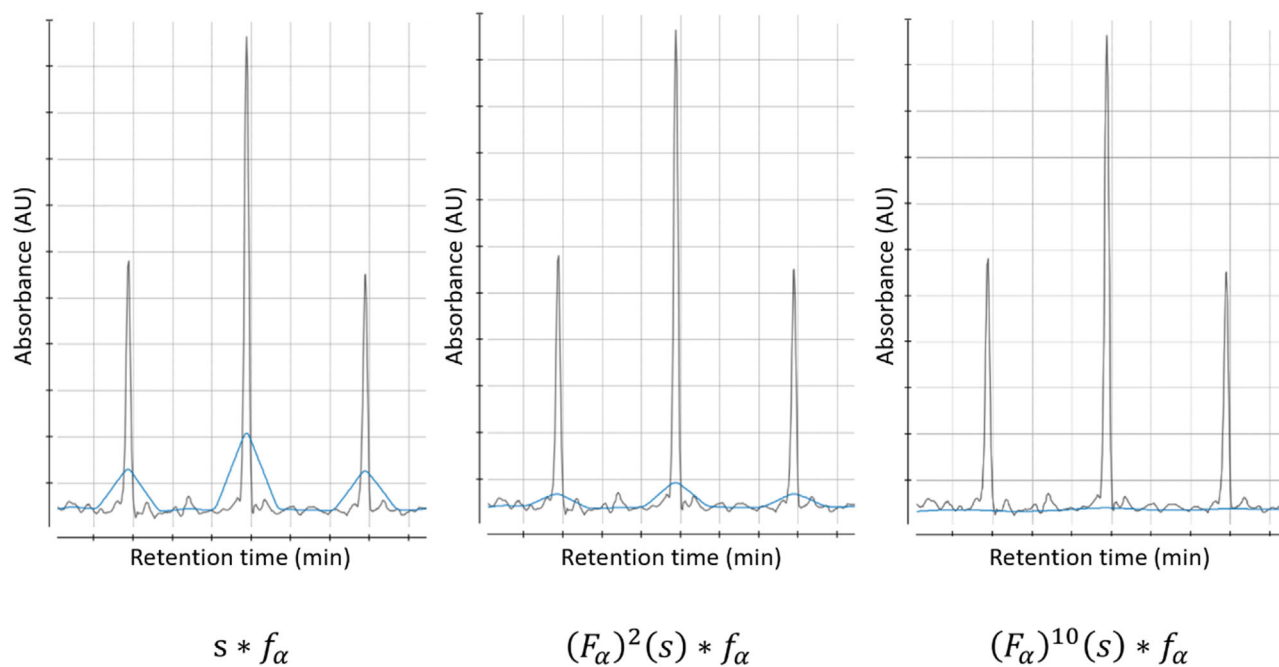


FIGURE 3 Influence of the number of iterations on the baseline estimation ($\alpha = 20$)